

**onyx** RoBa-Edition Extracting Forceps



**onyx**  
DESIGNED


**ZEPF**  
**roba**  
EDITION  
PATENTED

With the patented RoBa-Edition **HELMUT ZEPF** introduces a new generation of extracting forceps. Deduced from the **ZEPF Xcision** Instruments the RoBa-Edition has been especially developed in consideration of easy and gentle extraction.


The patented RoBa-Edition according to Dentist Beck is the consequential advancement of conventional extracting forceps with the advantages of the tapered deep-grip extracting forceps.

The modified beaks according to Dentist Beck fit exactly on the teeth which ensures a maximum grip in the appliance. These new beaks are available for all figures in upper and lower jaw (incisors, premolars, molars and wisdom teeth).


Due to the fact that all teeth show a convex crown contour (upper jaw: labial, buccal, palatal and in lower jaw: labial, buccal, lingual), the beaks have been developed under this anatomical actuality. The handle is a protected design from **HELMUT ZEPF**, which was developed in cooperation with Dr. Maty, Germany.

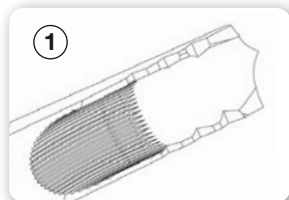


The bilow-grind claws radially into the teeth and avoids a "riding" in-between the tooth crown and the inner contour of the beak.

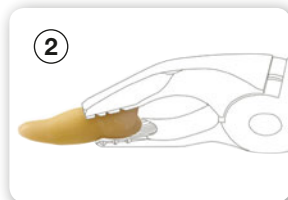

Concavely elaborated inner contours of the forceps beak fit in the convex tooth contour. With deep grip in the alveolus or on the crown, the RoBa Edition ensures a parallel and maximum grip in any situation. Root fracture almost can be excluded.


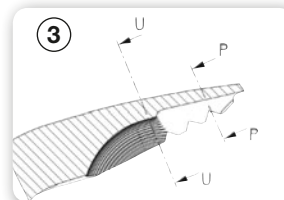


The different deeply elaborated inner contours ensure maximal adaptation on the teeth in different actualities. No tilting of the teeth in rotary and /or lifting movement.

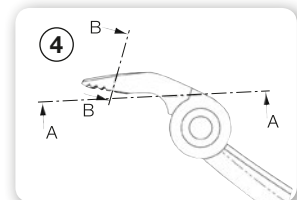

Tapered outside contour of the beak affords deep grip even subgingival.

scale 1:1

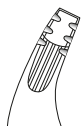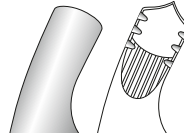						
						
3.2.1.   1.2.3.	5.4.   4.5.	5.4.   4.5.	7.6.	6.7.	8.   8.	8.   8.
<b>12.234.08Z</b>	<b>12.235.07Z</b>	<b>12.235.08Z</b>	<b>12.217.00Z</b>	<b>12.218.00Z</b>	<b>12.267.01Z</b>	<b>12.279.90Z</b>
<b>12.234.08ZD</b>	<b>12.235.07ZD</b>	<b>12.235.08ZD</b>	<b>12.217.00ZD</b>	<b>12.218.00ZD</b>	<b>12.267.01ZD</b>	<b>12.279.90ZD</b>
<b>12.234.08ZTI</b>	<b>12.235.07ZTI</b>	<b>12.235.08ZTI</b>	<b>12.217.00ZTI</b>	<b>12.218.00ZTI</b>	<b>12.267.01ZTI</b>	<b>12.279.90ZTI</b>



### Black Finish Coating

Design meets functionality. Users with highest demands will appreciate the elegant, black finish providing the instrument with a non-reflecting, extremely smooth and scratch-resistant surface. The article number is complemented by **TI**.



12.234.07ZTI (scale 1:1) PATENTED



3.2.1. | 1.2.3.

12.234.07Z

12.234.07ZD

12.234.07ZTI



### Diamond Version








All RoBa-Edition extracting forceps in the classic satin-metallic finished surface are available with diamond coating for a better grip.

The lower part of the ZEPF X vision Handle is gold-plated.

For ordering etc., just add to the item number the letter **D**.



12.234.07ZD (scale 1:1) PATENTED

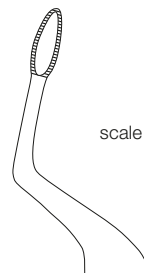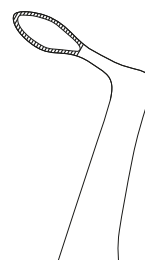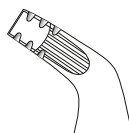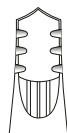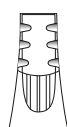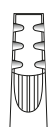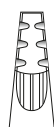

scale 1:1



Rescue-Line

Rescue-Line

W

W

3.2.1. | 1.2.3.

3.2.1. | 1.2.3.

5.4. | 4.5.

7.6. | 6.7.

7.6. | 6.7.

▶ ◀ ▶ ◀

▶ ◀ ▶ ◀

12.236.07Z

12.236.08Z

12.213.00Z

12.222.00Z

12.221.90Z

12.045.15ZS

12.051.15ZS

12.236.07ZD

12.236.08ZD

12.213.00ZD

12.222.00ZD

12.221.90ZD

12.045.15ZSD

12.051.15ZSD

12.236.07ZTI

12.236.08ZTI

12.213.00ZTI

12.222.00ZTI

12.221.90ZTI

12.045.15ZSTI

12.051.15ZSTI

**ZEPF**  
**roba**  
EDITION  
PATENTED

Dedicated to the long-standing employees Mr. Horst Roos and Mr. Reinhold Bacher, who worked for Helmut Zepf for more than 50 years since their time as trainee.

The **ZEPF** RoBa-Edition completes the innovative advancement of the famous **ZEPF** Extracting Forceps.

The RoBa Edition is developed under the aspects of a secure and gentle appliance in direct collaboration with practitioners of incomparable experiences.



AESTHETIC IS THE RESULT