

LANGER Universal Curette # L 17/18

in **ZEPF** *bionik* Handle



24.207.17L



LANGER Universal Curette # L 17/18

in **ZEPF BIONIK** Handle


The new Langer Curette # L 17/18 has the same bending as the Gracey # GRA 17/18 and thus is the perfect complement to our range of Universal Curettes.

Unlike the Gracey Curettes, Universal Curettes have two working ends / cutting edges, so that concretions on the tooth neck as well as inflamed tissue in the pockets can be removed simultaneously in one step.

The triple bend of the working ends provides optimal access, especially in the distal premolar and molar area and even when mouth opening is restricted.

The **ZEPF BIONIK** handle is designed according to the ergonomic needs of a curettage.

QUICKFIX

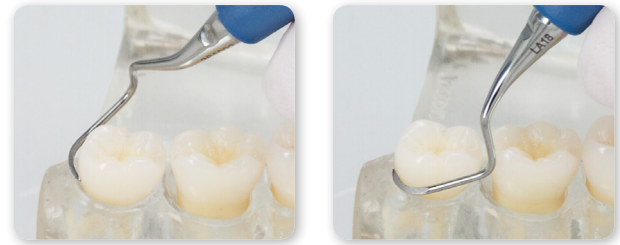

8-4 | 4-8
8-4 | 4-8

24.207.17L

Langer # L 17/18, Universal Curette for use on upper and lower molars and premolars, in **ZEPF BIONIK** handle cobalt blue, inserts exchangeable

Option:

All teeth from 1-8 in the upper and lower jaws are covered when the universal curette Langer # L 5/6 for the anterior region is used together with the new Langer # L 17/18 (24.207.17L) for premolars and molars.



3-1 | 1-3
3-1 | 1-3

24.201.05L

Langer, # L 5/6, Universal Curette for use on upper and lower front teeth, in **ZEPF BIONIK** handle yellow, inserts exchangeable

BIONIK ZEPF



reddot design award
winner 2010

BIONIK Universal Handle made of PEEK high-tech plastic material – guarantees an ideal power transmission with formerly unknown sensitivity. The handle is available in 10 different, fresh basic colors.

The exchangeable working tips inserted in the ergonomic **BIONIK** handle offer highest economy and best tactile handling.

AESTHETIC IS THE RESULT