

**onyx** Hand Instruments  
DESIGN



Based on the motto "WORK SMART – NOT HARD", selected instruments for dental surgery are now provided with **ZEPF onyx** Coating.

The black coated **onyx** Instruments, in the ergonomic **bionik** handle each, dispose of a durable, very hard and impact-resistant surface with excellent sliding qualities.

The extremely hard surface guarantees a long stability of the cutting edges, and a regular re-sharpening is no longer required. The black, scratch-proof surface is easy to clean and avoids unpleasant light reflections. Furthermore, a drying of impurities is avoided due to the smooth surface.



Raspatory round / flat



Raspatory, bent 90°, slim version. Very well suited for the first dissection of sinus mucosa, before mobilization with the blunt elevator



Molt Curette Combination of CM4 / CM7, # 4 Ø 3.5 mm / # 7 Ø 7.0 mm



Lingual Raspatory Combination Back Action for optimal lingual and palatal access



Tunnelling Instrument, excellent contrast in the soft tissues when tunnelling in the upper jaw

**24.748.52OX**

**24.747.03OX**

**24.204.35OX**

**24.747.07OX**

**46.040.01OX**



The **ONYX<sup>DESIGN</sup>** coating is applied to the carrier material in a more gentle way than conventional procedures, as there is less influence of heat. This coating is more durable, harder and impact-resistant. It also presents excellent sliding qualities favoring the care of the instruments. An incrustation of proteins etc. is not possible. Dirt can be wiped off. Even if the cutting edge does not appear black, the coating is nevertheless integrated in the surface.

**24.751.522MOX**



Tunneling Instrument, excellent contrast in the soft tissues when tunneling in the lower jaw

**24.751.520WOX**



Weng Papillex – the papilla elevator with the special blend, to prepare the papilla in a smart way

**26.182.320X**



Periotome Combination Instrument, optimal surface hardness and best contrast\*

**24.751.178OX**



Back -Action Raspatory as bone scraper with dissection tip. This tip serves to present individual layers of tissue in anatraumatic way (dissection) during the preparation of the soft tissue

**41.751.011LUOX**



Lucas Scraper, delicate version, to clean the extraction alveolus and to remove granulation tissue

**24.751.523MOX**



Tunneling Instrument, excellent contrast in the soft tissues when tunneling in the lower jaw

**24.751.521WOX**



Weng Papillex – the papilla elevator with the special blend, to prepare the papilla in a smart way

**26.182.330X**



Periotome Combination Instrument, optimal surface hardness and best contrast\*

**24.751.512OX**



Back -Action Raspatory as bone scraper with dissection tip. This tip serves to present individual layers of tissue in anatraumatic way (dissection) during the preparation of the soft tissue

**41.751.012LUOX**




Lucas Scraper, delicate version, to clean the extraction alveolus and to remove granulation tissue



## Application Pictures

from Dr. Claus Fuchs, Konstanz-Litzelstetten / Germany

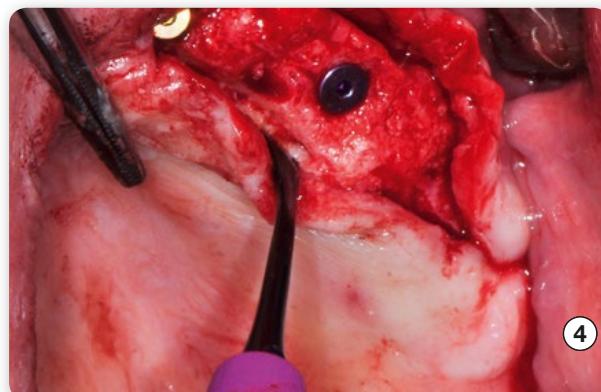

AESTHETIC IS THE RESULT